


SKU: 820534 UPC: 817475001001

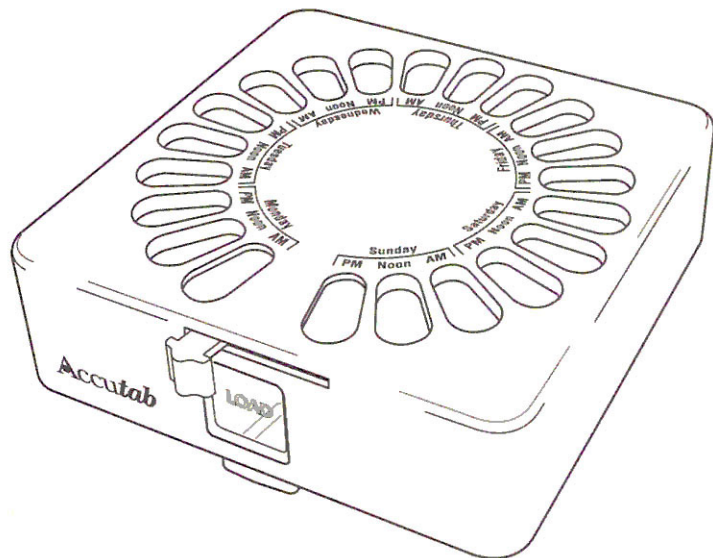

Accutab is NOT childproof.

Congratulations on your purchase of Accutab. Accutab is an innovation in dosing for prescription medications and vitamin supplements. It simplifies the organization of weekly and multiple medications, and helps ensure correct dosing.

Accutab simplifies the task of medication management for persons of all ages. Accutab has an added bonus for arthritis patients, whose condition makes it painful to open medicine bottles multiple times each day. The concept for the Accutab dispenser was developed by an arthritis patient.

Accutab can be used with an electronic medication reminder such as the e-pill Multi Alarm HOME (SKU 759023), e-pill CADEX Alarm Watch (SKU 952431 or 952433), or any Smartwatch.

Visit www.epill.com for more information.

How to Mount your Accutab™ Using Included Screws

Permanently mounting your e-pill weekly medication dispenser before its first use is recommended. Your new e-pill dispenser comes with an easy to install bracket (screws are included) designed to mount inside a kitchen cabinet (the recommended child-proof location) or under a shelf.


A screwdriver is the only tool required for installation. A pencil is useful for marking the mounting location. The bracket is attached to the mounting surface with four screws provided.

1. Find a convenient and suitable location to mount your Accutab. Some oral tablets and capsules react adversely to high humidity. Mounting directly over sinks or in bathrooms is not recommended. **Accutab is not child proof** and we recommend mounting inside a cabinet, out of reach of children.

2. Insert Accutab into its bracket and hold the assembly in place, making sure there is clearance for removing

3. While holding the assembly in place, mark the screw locations on the mounting surface through the holes in the bracket with a pencil or pen.


4. Remove Accutab from the bracket and attach the bracket to the mounting surface with the four screws provided. Once the bracket is attached, insert Accutab and test its operation.

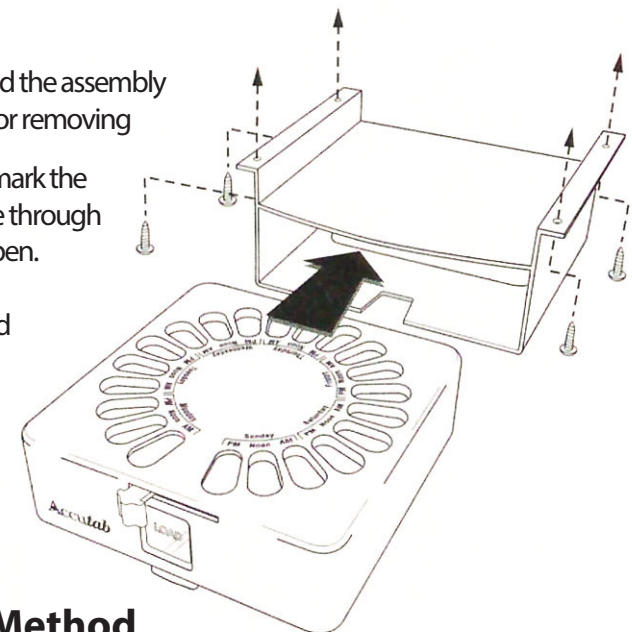

How to Mount your Accutab™ - Alternative Method

Alternatively, you can use 3M™ VHB™ Double-Sided Tape to mount your Accutab and avoid the use of a drill. **Follow Steps 1 through 3 above.** Instead of drilling holes, use the pencil marks as guides.

Adhere the VHB™ Tape inbetween the holes, facing the surface. Peel off the protective backing and carefully line up the holes and pencil marks, then simply push the unit up against the surface to mount your bracket. Apply even pressure for about 30 seconds.


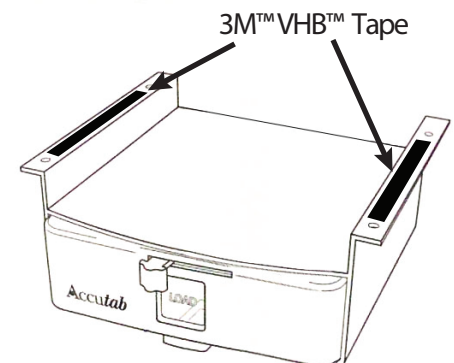


Loading and Operation of Your Accutab™

Accutab is intended to provide up to three doses per day starting with the Monday AM dose and continuing through Sunday PM. Accutab should be loaded after the Sunday PM dose and before the Monday AM dose. Accutab should only be loaded when the word LOAD appears in the Dose Window (which shows the dose last taken).

1. Ensure the word LOAD is visible in the dose window. If LOAD is not visible move the red dispensing lever from left to right until LOAD is visible.
2. When loading, remove Accutab from its bracket and place it on flat surface.
3. Load your medications in the appropriate compartments according to your dosing schedule. For each Day of the Week, three compartments are provided and are labeled AM, NOON, and PM. You may fit up to 25 aspirin-sized tablets in each dosage compartment. **See the next page for specific filling instructions based on how many times you take medication.**
4. After loading, insert Accutab into its bracket. Accutab is now ready for operation. **The dose that appears in the clear window is the dose LAST TAKEN.**

Operating the Accutab is very simple. Place your hand or a container ("medicine cup") under the dispensing opening located on the bottom of Accutab. Moving the red dispensing lever from left to right one time will dispense one dose at a time. The dose will fall from the dispenser opening. The compartment day/time label viewed through the dose window is a convenient reminder of the LAST DOSE TAKEN (last dose dispensed). The clear compartment walls allow you to see that all medications have fallen from the compartment, and that no medications are left in the tray.

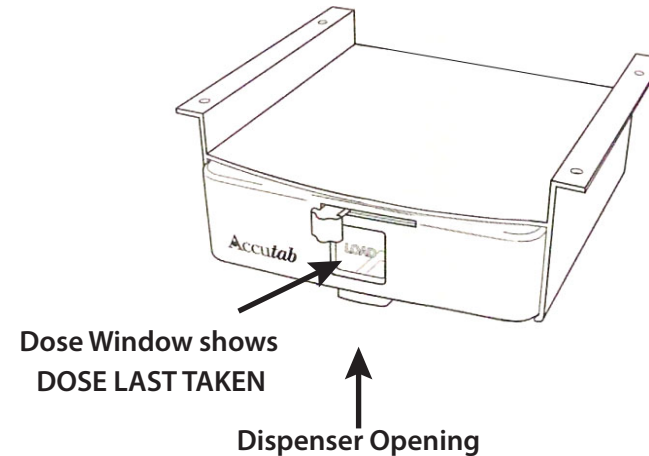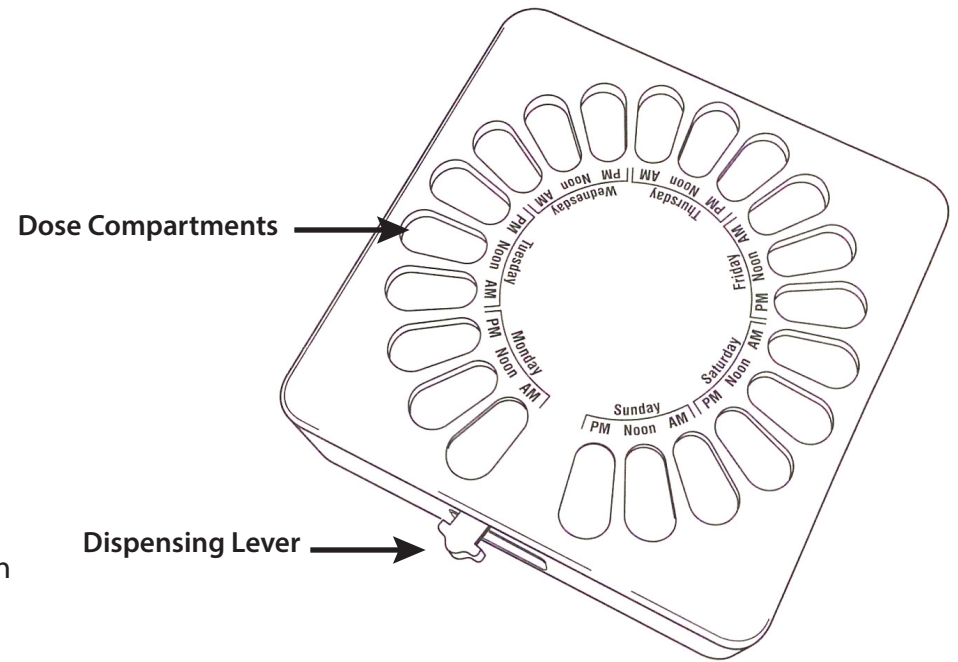


**Accutab is
NOT childproof.**

One Year Limited Warranty

e-pill, LLC (e-pill Medication Reminders) warrants this product to be free from defects in material and workmanship for a period of one year from the date of purchase. If this product is found to be defective, it will be replaced at no cost to you. You may contact e-pill, LLC in the US and Canada toll free by calling 1-800-549-0095 or +1-781-239-2941. This warranty does not apply to normal wear or to damage resulting from abnormal use, misuse, abuse, neglect, or accident. e-pill, LLC will not be liable for any accidental, consequential or special damages, however caused, for breach of warranty or otherwise.

This e-pill Medication Reminder is only a supplement to any memory system that you currently use to remind yourself to take medications. The user maintains full responsibility to administer medication at the appropriate times.


e-pill, LLC assumes no responsibility for any incidental or consequential damages, including but without limitation, damages or injuries resulting from malfunction of the product. e-pill, LLC's liability shall in no circumstances exceed the amounts paid by the user for the e-pill Accutab.

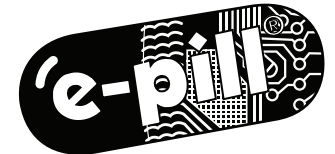

e-pill® | Medication Reminders
49 Walnut St, Bldg 4
Wellesley, MA 02481, USA
1.800.549.0095 | epill.com

How to Fill the Accutab™ if you take Medications *Once Per Day*

1x

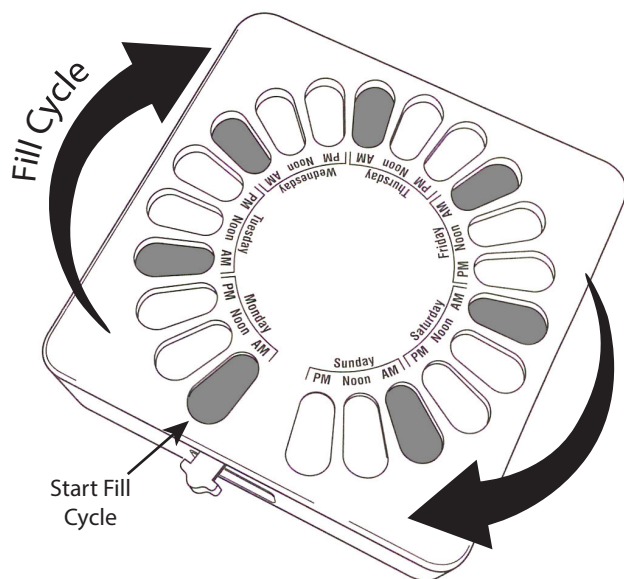


■ Fill these compartments

The user will have to move the dispensing lever 3 times to get to the daily dose.

How to Fill the Accutab™ if you take Medications *Twice Per Day*

2x

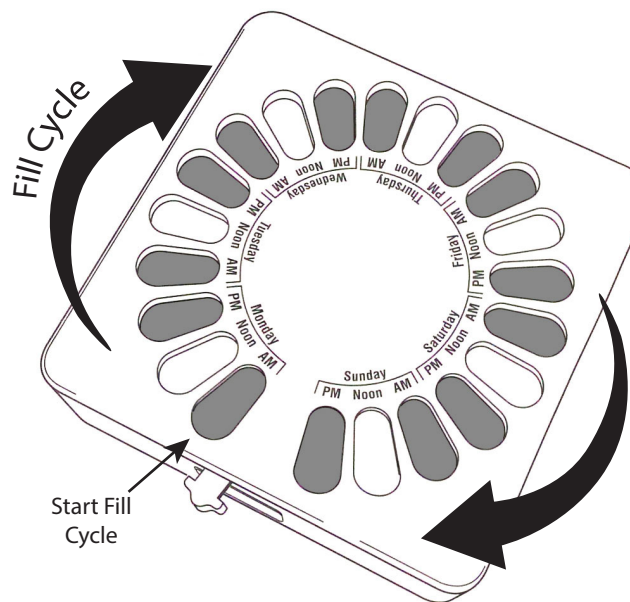


■ Fill these compartments

The user will have to move the dispensing lever 1 time to get to the first daily dose.

They will have to move the dispensing lever 2 times to get to the second daily dose.

How to Fill the Accutab™ if you take Medications *Three Times Per Day*

3x

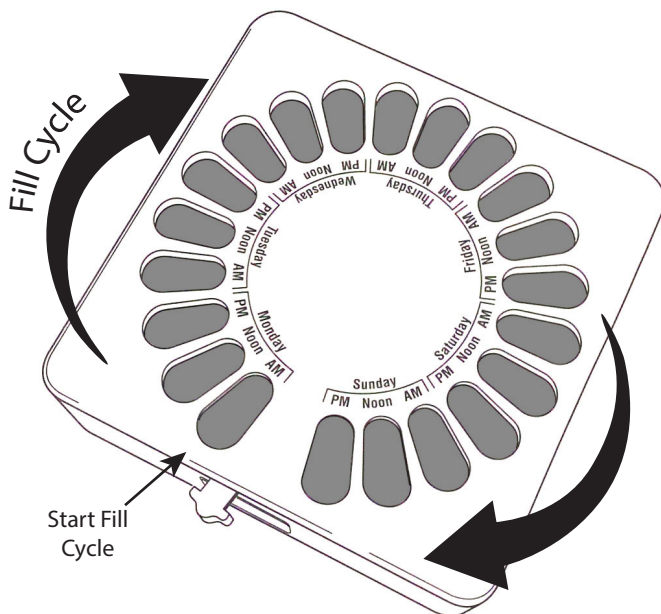

■ Fill these compartments

The user will have to move the dispensing lever 1 time to get to each daily dose.

The LAST TAKEN Indicator



The unique LAST TAKEN Indicator (dose window) shows the patient or caregiver when medications were last taken. It also allows the user to see that all medications have fallen from the compartment and that there are no medications left in the tray.

For the User of the e-pill® Accutab™

What to Do When it is Time to Take Your Medication

- 1)** Place your hand or a container ("medicine cup") under the dispensing opening located on the bottom of Accutab.
- 2)** Move the **red** dispensing lever from **left** to **right** one time to dispense one dose at a time. The dose will fall from the dispenser opening.
- 3)** Take the medication that has fallen into the the palm of your hand or medicine cup.

That's it! Simple to use for everyone!

The compartment day/time label viewed through the dose window is a convenient reminder of the **LAST DOSE TAKEN** (last dose dispensed). The clear compartment walls allow you to see that all medications have fallen from the compartment, and that **no medications are left in the tray.**

The Accutab holds up to 25 aspirin-sized pills per dose.

Components of the e-pill® Accutab™

- 1)** The Accutab Weekly Organizer
- 2)** The Mounting Bracket
- 3)** Mounting Screws
- 4)** Mounting Tape
- 3)** A Daily Pill Box for On- The- Go Use - If you are away and your medications need to be taken, use this portable pillbox.

