


MEDICATION MANAGER

Easy-to-use pill dispenser with clock & talking reminder

Instruction Manual

INTRODUCTION


Congratulations on your purchase of a Timex Medication Manager-Model number 80041 for the Daily Medication Manager or Model number 80040 for the Weekly Medication Manager. The Medication Managers are designed to help you take your medications on a timely manner with ease and confidence. Please read the instruction manual in full before using your unit. If you have any questions, you may contact us at 1-800-422-0360.

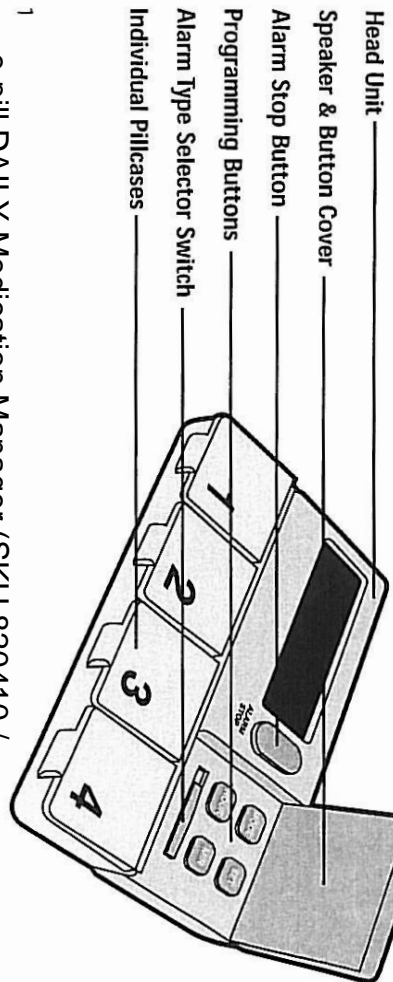

- 1 e-pill DAILY Medication Manager (SKU 820410 / www.epill.com/medmanager.html)
- e-pill WEEKLY Medication Manager (SKU 820470 / www.epill.com/weeklymed.html)

SETTING THE CLOCK

Slide pill box to left. Lift the speaker & button cover on the face of the timer and press the **MODE** button to choose what function you want to set. While in the Time mode, press and hold the **SET** button for three (3) seconds.

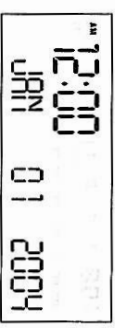

AM 12:00 will flash.

Press **HOUR** button to increase hour. Once desired time is reached press **MIN** button to increase minutes. Once desired time is reached, press **SET**. The time will stop flashing and is now set.

The month abbreviation will now flash. To set the correct month, press the **HOUR** button until correct month shows then press **SET**.

The day of the month will begin to flash. Press the **HOUR** until the correct date shows the press **SET** button. The year will begin to flash. Press **HOUR** until the correct year shows then press **SET**.


Note: If you pass desired time/date/month you will need to keep going forward by hitting the **HOUR** button until the correct information is displayed. Press **SET** and continue.

Note: If you are making any changes to the clock and 30 seconds lapse without any buttons being pushed, it will return to normal clock mode.



INSERTING BATTERIES

To insert batteries, slide open battery door and insert two (2) AAA size batteries making sure the negative (-) and positive (+) terminals are placed correctly. The product comes with 2 LR44 batteries (see illustration). To activate, pull and remove the insulation tab. The LR44 batteries are to hold memory when AAA's are running low or are being replaced. They are replaceable.

A battery indicator will show on the display when you need to replace batteries.

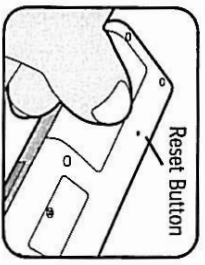

2 AAA Battery Replacement

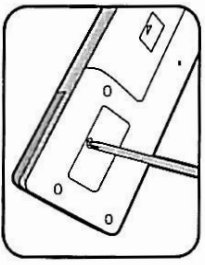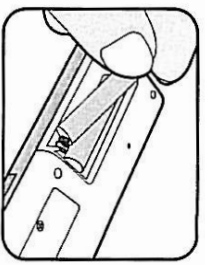

2 LR44 Battery Replacement

SETTING THE TIMER

You may now choose up to 4 time settings – one for each pillbox compartment. You will also have the option to set 2 additional alarms as personal reminders. These may be used to remind you to take a blood pressure reading, blood glucose test, etc...

Press **MODE** key.



ALARM OFF will appear (see illustration). This will be the reminder time for pillbox number 1 and the number will be flashing.

Press and hold **SET** key for three (3) seconds. --- 1 will flash. Press **HOUR** button to set the correct hour, then press **MIN** to set the correct minute. Once desired time shows, press the **SET** button. Time will stop flashing and **ALARM ON** will appear.

Press **MODE** key again to set the alarm time for pillbox number 2. Repeat process above to set alarms 2 through 4.

To change a dosage time, simply select the dosage that needs changed and repeat the above set process.

Note: If only 1, 2 or 3 reminders are needed, the boxes should be set in sequential order. The timer will not allow you to set 1, 2 and 4 box. You will see **ALARM SET ERROR**.

4

SETTING ALARM TYPE

There are three (3) alarm types to choose from:

Beeping alarm – Indicated by a bell.



Voice alarm – "TAKE PILLBOX 1" indicated by "Voice" and a bell.

Silent – No alarms heard, only prompted by display.

The alarm type can be changed at any time. Simply slide the button to the far left for the Beeping alarm, Center position for Voice alarm or move to the far right for the visual or silent mode.

MISSED PILL INDICATOR

Once all the settings are in place you will be reminded when it is time to take your dosage.

You will either hear the voice, beeping or see the visual reminder and the pillbox that you are to take will flash and have a box around the number. It will also scroll across the display "Take Pillbox 1".

Push the **ALARM STOP** button next to the display window to indicate that you have taken your medication. If the **ALARM STOP** is not pressed, you will be reminded once every minute for 3 minutes. If dosage has still not

6

When you need to delete a dosage time or make several changes for your dosages, you will need to make those changes in reverse order of setting, start with box 4, then 3, then 2 and finally 1. You can also just press the small **RESET** button on the back of the timer using a small tip pen or paper clip end (see battery changing illustrations). It will clear all set information and allow you to start over.

2 additional timers

After completing and setting the alarm times for each of the 4 compartments, you will have the option of 2 additional timers.

Press **MODE**. Display will be what is illustrated at right.


Press and hold **SET** key for three (3) seconds. --- will flash.

Press **HOUR** button to set the correct hour, then press **MIN** to set the correct minutes. Once desired time shows, press the **SET** button. "Blood Pressure Alarm On" will scroll across the screen, stopping on Alarm On.

To set BP2, press **MIN** and repeat steps above to set optional BP2 alarm.

Note: When setting any of the alarms, the minimum time between must be at least 10 minutes. An **ALARM SET ERROR** will scroll across the screen if the minimum is not observed. Press and hold **SET** and repeat the setting process as listed above.

5

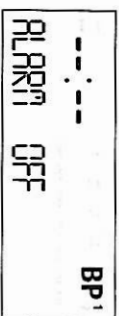

been taken and **ALARM STOP** not pressed, the dosage will be considered "Missed". That dosage number will then be shown with an X (refer to illustration) and the Message "MISSED PILLBOX 1" will display on screen until it is time for next reminder.



HISTORY

After a minimum of 1 full day of reminders, you will be able to view **HISTORY** and can see what doses were missed. One full day is measured by all dosages before 11:59PM.

Press **MODE** key 7 times. **HISTORY** will be shown on screen. By pressing **SET** you will be able to see last 7 days history one day at a time. You will see all missed and taken dosages during that day.




7

TRAY INSTALLATION

If you have purchased the Weekly Medication Manager Model #80040, you may want to mount the tray to a vertical surface such as a wall or cabinet door. The tray may also be used on a desk or countertop.

The tray comes with 2 screws for the installation. Determine the desired location and mark where the top screw will go. Install the top screw but do not tighten fully. Now mark where the lower screw will go. Remove the tray and install the lower screw. Reinstall the tray and tighten the screws.

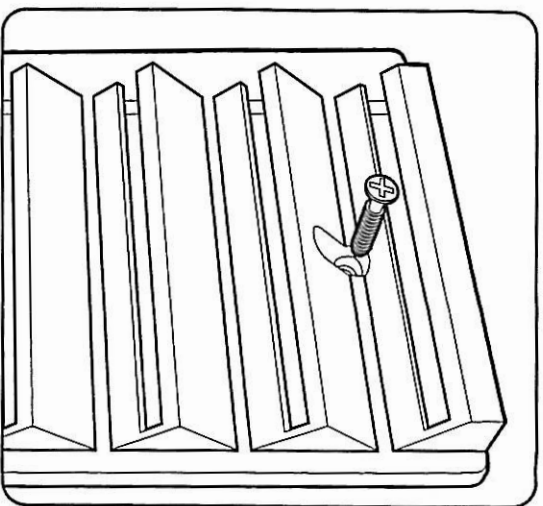

8

Limited Warranty

MEDport takes pride in the quality and workmanship crafted into every product we make. Accordingly, MEDport will repair or replace at our option any product which is defective in material or workmanship for a period of one year from the date of purchase. This warranty does not cover damage caused by accident, negligence or improper care. Product should be returned to MEDport with proof of purchase and postage prepaid along with an explanation of the problem.

TIMEX is a trademark of Timex Corporation.

Manufactured and Distributed under license by:

MEDport, LLC, 23 Acorn Street, Providence, RI 02903

Customer Service 1-800-422-0360

8:00 AM - 8:00 PM EST -Monday through Friday

www.timexhealthcare.com

Email: info@medportllc.com

Made in China

Detailed user information about our Medication Managers can be found at www.timexhealthcare.com.

9

To re-order please visit or call:

www.epill.com 1-800-549-0095

e-pill DAILY Medication Manager (SKU 820410 / www.epill.com/medmanager.html)

e-pill WEEKLY Medication Manager (SKU 820470 / www.epill.com/weeklymed.html)